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제품구성 :
509-2 One-Fil Putty 0.25g x 4 Syringes

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< 3.5h icon"/>

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Australian Endodontic Journal
The Journal of the Australian Society of Endodontology Inc,
the Australian and New Zealand Academy of Endodontists and
the Asian Pacific Endodontic Confederation


Aust Endod J 2021

ORIGINAL RESEARCH

Physicochemical properties and cytocompatibility of newly developed calcium silicate-based sealers

Min-Gyu Park, DDS¹; In-Ryoung Kim, PhD^{2,3}; Hyung Joon Kim, PhD^{4,5}; Sang Won Kwak, DDS, MS, PhD^{1,2}; and Hyeon-Cheol Kim, DDS, MS, PhD^{1,3}

Cell Attachment Evaluation (One-Fil)

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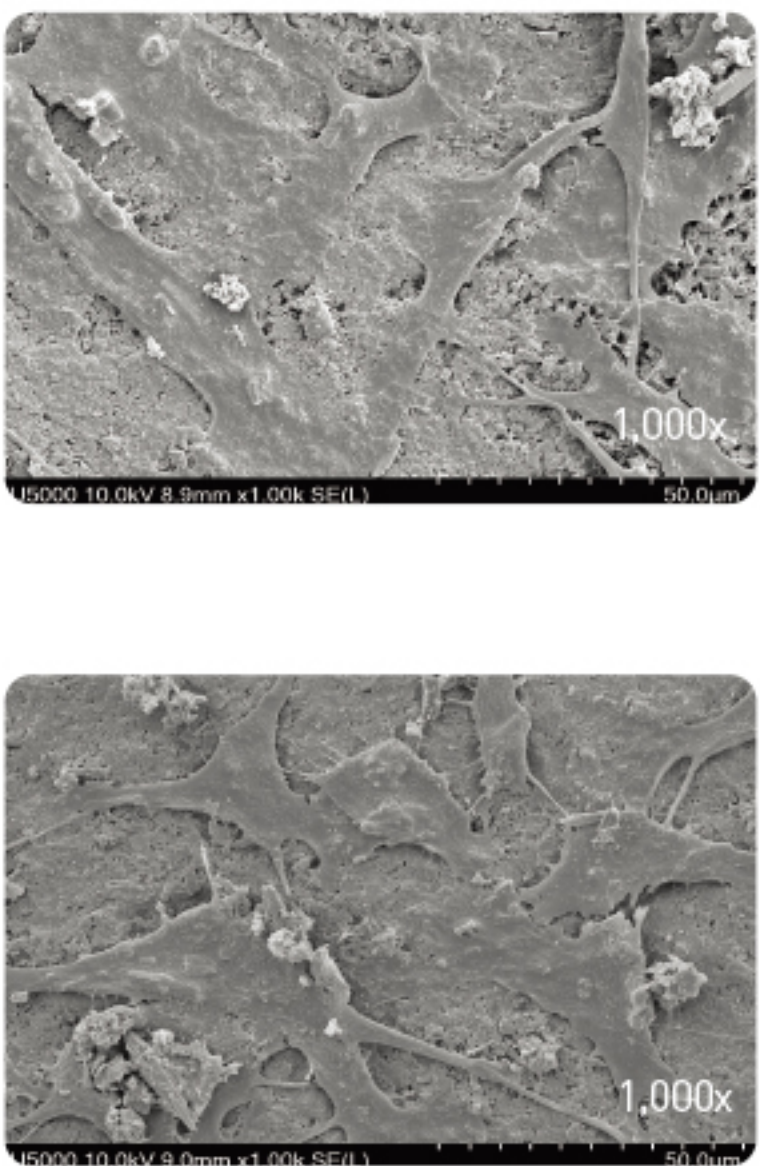BASIC RESEARCH - TECHNOLOGY

Chemomechanical Properties and Biocompatibility of Various Premixed Putty-type Bioactive Ceramic Cements

Minju Song, DDS, MS, PhD,¹ So-Mang Lee, DDS, MS,¹ Ji-Young Bang, BS, MS,¹ Ruben H. Kim, DDS, PhD,² Sang Won Kwak, DDS, MS, PhD,¹ and Hyeon-Cheol Kim, DDS, MS, PhD¹

ABSTRACT
Introduction: This study aimed to evaluate the chemomechanical properties and biocompatibility of recently introduced premixed putty-type bioactive ceramic cements (PPBCs). Methods: Including ProRoot MTA (PMTA) as a control, BC RRM fast-set putty (BCPT), Well-Fix PT (WFPT), One-Fil PT (OFPT), and Endocem MTA premixed (ECPM) were compared to evaluate setting time, radiopacity, pH change, and microhardness. Biocompatibility on human dental pulp cells was compared using CCK-8 assay. Mineralization potential was evaluated using alkaline phosphatase activity, Alizarin Red S (ARS) staining, and quantitative real-time polymerase chain reaction with odontogenic gene marker. For data analysis, 1-way analysis of variance and Tukey's post hoc test were used at the significance level of 50%. Results: Among the PPBCs, BCPT presented the longest (52 ± 27) setting time (minutes) and others showed significantly shorter time than PMTA (34 ± 22) (P < .05). WFPT (6.20 ± 0.54) and OFPT (6.82 ± 0.53) showed significantly higher radiopacity values (mmAl) and others showed similar value compared with PMTA (P > .05). All PPBCs showed high alkaline pH from fresh materials and tended to increase according to time periods from 30 minutes to 12 hours. ECPM showed the highest value of microhardness (81.62 ± 5.93). WFPT showed similar, and others showed lower than PMTA (P < .05). All PPBCs showed biocompatibility in CCK-8 assay. All PPBCs showed similar or better value compared with PMTA in ALP and ARS staining, and ALP and DSPP marker expression (P < .05). Conclusions: The PPBCs showed clinically acceptable chemomechanical properties and favorable mineralization potential. (J Endod 2023;49:1713-1721.)

KEY WORDS
Biocompatibility; mineralization potential; microhardness; premixed putty-type bio ceramics;

SIGNIFICANCE

The recently introduced premixed putty-type bio ceramics, which have the advantages of uniform consistency and user-friendly manipulation, showed clinically acceptable chemomechanical properties and favorable mineralization potential.

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