

Any-Com[®] Seal

Effective Pit & Fissure Sealant Containing 10-MDP and MBG

Powerful Protection for Pits and Fissures



Same Procedure, Enhanced Results

Enhanced bonding through functional monomers
Enamel remineralization promoted by MBG
(Mesoporous Bioactive Glass)

A Reliable Tooth Protection System

37% phosphoric acid gels available in different viscosities:
Dentto-Etch Plus, Any-Etch
Desens Bio-Care gel for additional protection in the clinic and at home

Any-Com® Seal

Fast and Comfortable Non-Invasive Fissure Sealing for Children

Any-Com® Seal is a light-cured composite resin-based sealant that provides long-term protection through continuous fluoride release. Formulated with 10-MDP and MBG (Mesoporous Bioactive Glass), the sealant enhances adhesion to tooth surfaces and stimulates hydroxyapatite (HA) formation, a major component of enamel.

As a result, the sealed surface and strengthened enamel effectively protect teeth from cariogenic bacteria.


Indications

- Sealing of pits and fissures

Benefits

- Enhanced adhesion and remineralization: 10-MDP, MBG
- High wear resistance and strength with 42% inorganic filler
- Smooth dispensing and precise application with optimal viscosity and flowability
- Excellent penetration even in complex occlusal anatomy and upper molars thanks to thixotropic behavior

Fluoride Release

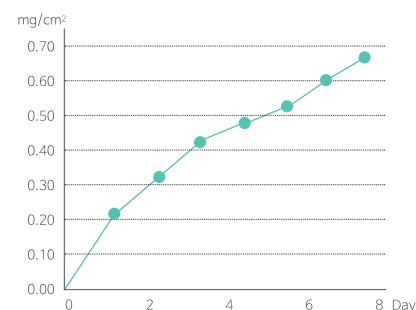

* Source: MEDICLUS R&D

Step by Step Process



Scientific Publication



Journal of Dentistry
Volume 160, September 2025, 105969

Enhancing dental sealant performance:
Effects of mesoporous bioactive glass and
10-MDP on adhesion and remineralization

Nara Yun ^{a,1}, Hyun-Jung Kim ^{a,1}, Jiyoung Kwon ^b, Sun-Young Kim ^c, Duck-Su Kim ^{a, R, ✉}

Delivery Forms



#308-1 Any-Com® Seal

Any-Com Seal 2g x 2ea, White

Flex Tip x 20ea

MEDICLUS

Headquarter | No.1210, 134, Gongdan-ro, Heungdeok-gu, Cheongju-si, Chungcheongbuk-do, 28576, Korea | +82 43 211 2877

Sales Office | No.206 Brown Stone, 464, Cheongpa-ro, Jung-gu, Seoul, 04510, Korea | +82 10 8385 2877

www.mediclus.co.kr | sales@mdclus.com

