

Any-Core®

FDA


DUAL CURED, CORE BUILD-UP COMPOSITE RESIN

Any-Core is dual cured composite material for use in the fabrication of core build-up and post cementation.


Advantages

- Optimum viscosity (flow) & flowability
- Excellent compressive strength (330MPa) and flexural strength (160MPa)
- Preparation similar to dentin
For A3 color, you can easily prepare the treated tooth, because the hardness of our core material is similar to that of dentin.
- light-curing time 20sec / 2mm
- Chemical curing time : 5min (intra oral)
- Low polymerization shrinkage
- Excellent radiopaque

Physical properties

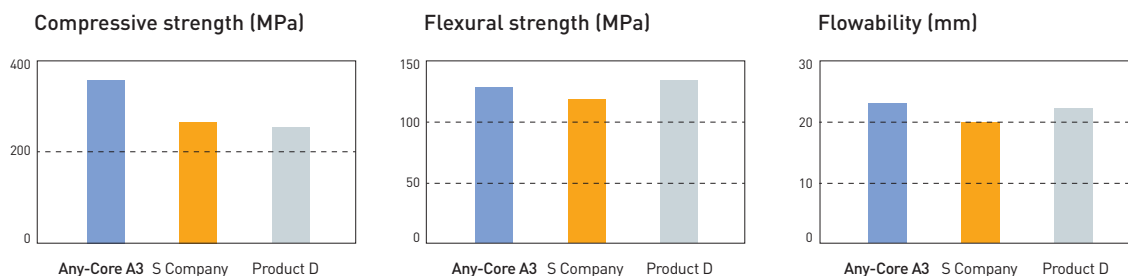

Application

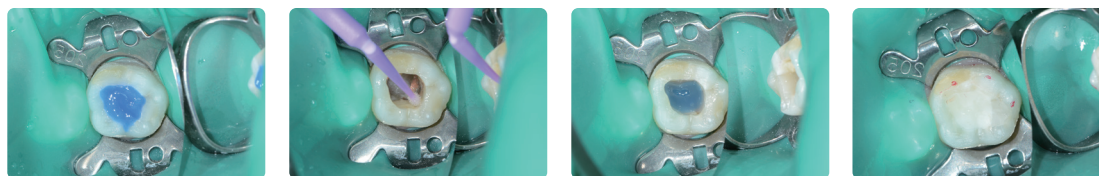


1. Etching

2. Bonding

3. Apply Any-Core Blue or A3

Results

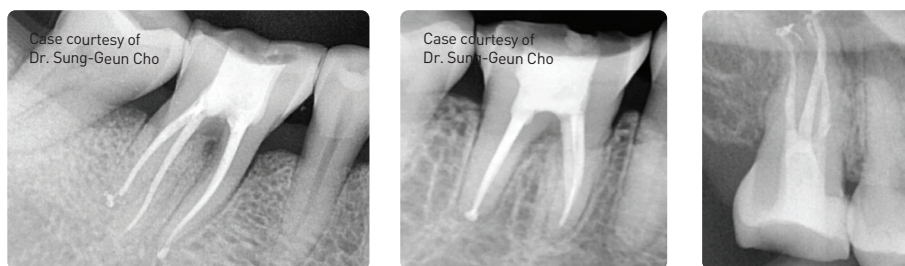

Case courtesy of Dr. Sung-Geun Cho

Case courtesy of Dr. Sung-Geun Cho



COMPOSITE RESIN FOR CORE BUILD-UP FDA

SHADE A3, White, Blue

CONTENTS

# 801-1 Any-Core Kit A3 9g x 2ea	# 802-2 Any-Core Refill White 9g x 1ea
# 801-2 Any-Core Refill A3 9g x 1ea	# 803-1 Any-Core Kit Blue 9g x 2ea
# 802-1 Any-Core Kit White 9g x 2ea	# 803-2 Any-Core Refill Blue 9g x 1ea